

Kineflex lumbar artificial disc versus Charité lumbar total disc replacement for the treatment of degenerative disc disease: A randomized non-inferiority trial with minimum of 2 years' follow-up

Kenneth Pettine and Andrew Hersh

Int J Spine Surg 2011, 5 (4) 108-113
doi: <https://doi.org/10.1016/j.esas.2011.07.003>
<https://www.ijssurgery.com/content/5/4/108>

This information is current as of November 2, 2024.

Email Alerts Receive free email-alerts when new articles cite this article. Sign up at:
<http://ijssurgery.com/alerts>

Kineflex lumbar artificial disc versus Charité lumbar total disc replacement for the treatment of degenerative disc disease: A randomized non-inferiority trial with minimum of 2 years' follow-up

Kenneth Pettine, MD ^{a,*}, Andrew Hersh, MD ^b

^a Loveland, CO

^b David Grant USAF Medical Centre, Travis AFB, CA

Abstract

Background: The Kineflex lumbar artificial disc replacement device (SpinalMotion, Mountain View, California) is a semiconstrained, posterior center of rotation, metal-on-metal intervertebral disc prosthesis. We performed a prospective, randomized, non-inferiority trial comparing the Kineflex Disc with the Food and Drug Administration (FDA)–approved Charité device (DePuy Spine, Raynham, Massachusetts). Our objective was to evaluate the Kineflex Disc's safety and efficacy using validated outcomes measures—the visual analog scale (VAS) and the Oswestry Disability Index (ODI).

Methods: Sixty-four patients were randomized to receive either the Kineflex Disc or Charité device and were then followed up for up to 3 years. Patients completed VAS and ODI questionnaires and were evaluated clinically and radiologically for complication or device failure. Results were analyzed in terms of change in mean VAS score and ODI from baseline, as well as with a comparison of clinical success as defined by FDA investigational device exemption criteria. Non-inferiority was defined as a difference of less than 18 points in the VAS score and difference of less than 10 units on the ODI scale, in keeping with a previously established minimum clinically important difference.

Results: The mean improvement for the Kineflex Disc group at 24 months was 56.80 for the VAS score and 37.30 for the ODI. Similarly, the mean improvement in the Charité group was 54.43 for the VAS score and 38.40 for the ODI. At 2 years of follow-up, no difference was found in VAS scores between the two groups. The Kineflex Disc group was therefore found to be non-inferior (mean difference, 2.37; 95% confidence interval, –12.5 to 17.3; $P = .004$). In addition, at 24 months, 83% of patients in the Kineflex Disc group and 85% of patients in the Charité group met FDA-defined criteria for clinical success, with no difference between groups ($P = .802$).

Conclusions: This level I evidence shows the Kineflex Disc to be non-inferior to the Charité device in terms of pain reduction (VAS score) and FDA-defined clinical success at 24 months' follow-up. Both devices showed a high degree of safety.

© 2011 SAS - The International Society for the Advancement of Spine Surgery. Published by Elsevier Inc. All rights reserved.

Keywords: Lumbar artificial disc replacement; Charite; Kineflex

Lumbar artificial disc replacement (ADR) devices have been in use in the United States for several years and have generated a great deal of controversy in both the scientific literature and lay literature.^{1–3} Total disc replacement has previously been shown to be non-inferior to spinal fusion in terms of pain and disability reduction.^{4–7} Replacement, however, offers the potential advantage of maintenance of normal motion and reduced forces on adjacent vertebral levels, which may result in

adjacent-level disease.^{8–11} Currently, 2 lumbar ADR devices are Food and Drug Administration (FDA) approved: the Charité prosthesis (DePuy Spine, Raynham, Massachusetts) and ProDisc-L prosthesis (Synthes, West Chester, Pennsylvania). Several newer devices, including the Kineflex Disc (SpinalMotion, Mountain View, California), are now seeking to gain FDA acceptance. This newer technology hopes to improve wear profiles and further decrease facet stress and shear forces while maintaining normal or near-normal range of motion.

We compared the safety and efficacy of the Kineflex Disc with the Charité artificial disc. The Charité device gained FDA approval in 2004 after the completion of several clinical trials, including a 2-year multicenter study

* Corresponding author: Kenneth Pettine, MD, 3810 North Grant Ave, Loveland, CO 80538; Tel: 970-669-8881; Fax: 970-669-4200.

E-mail address: kpettine@spinerevolution.com

Table 1
Inclusion criteria

Age between 18 and 60 years
 Evidence of degenerative disc disease at either L4–5 or L5-S1 (only 1 level) with radiographic evidence (CT, MRI, plain film, flexion/extension films, myelography, discography, and so on) of mild to moderate osteophyte formation of vertebral endplates, loss of disc height >2 mm when compared with adjacent level, herniated nucleus pulposus, loss of water content on MRI (black disc on T₂-weighted image), or vacuum phenomenon
 History of back and/or radicular pain that is severe, ongoing, and recurrent
 At least 6 months of prior conservative therapy for discogenic back pain and/or prior nucleolysis, nucleoplasty, discectomy, or laminotomy (without accompanying facetotomy) at study treatment level
 ODI showing moderate disability
 VAS showing moderate disability

Abbreviations: CT, computed tomography; MRI, magnetic resonance imaging.

where it was shown to be non-inferior to anterior lateral interbody fusion (ALIF) with BAK cages and iliac crest autogenous graft.^{12–15} Longer-term follow-up data subsequently published confirmed that non-inferiority and preservation of motion were maintained over time, and a reanalysis of the original data in 2007 suggested possible superiority throughout the 2-year follow-up.^{16–18}

The Kineflex Disc is a 3-piece chrome-cobalt-molybdenum device, with a semiconstrained core, posterior center of rotation design, and mobile center of rotation. Its design allows 12 degrees of movement into flexion, extension, and left-/right-side bending. A retaining ring on the inferior endplate limits excursion of the captured sliding core to 2 mm in all directions and prevents dislodgment of the sliding core. The disc is introduced as a single unit with fin fixation to the vertebral endplate allowing for increased ease of implantation and improved reliability in positioning.¹⁹ A non-inferiority design was used because these features are notable only if equivalent safety and efficacy of the Kineflex Disc device can be established.

Methods

This study is a prospective, randomized, controlled non-inferiority trial assessing the safety and efficacy of the Kineflex Disc to relieve symptoms associated with degenerative disc disease at either the L4–5 or L5-S1 levels. Institutional review board approval was secured before initiation of the study. Informed consent was obtained from all patients, who were then randomized to receive either a Kineflex Disc (n = 33) or Charité disc (n = 31). All patients were blinded to their device assignment until after surgery. Each was evaluated preoperatively and postoperatively, as well as at 6 weeks and 3, 6, 12, 24, and 36 months. Patients completed visual analog scale (VAS) and Oswestry Disability Index (ODI) questionnaires and were monitored clinically and radiologically for device failure, complications, and need for reoperation. Inclusion and exclusion criteria were similar to those used in the original ADR-versus-fusion trials (a full review of inclusion and exclusion criteria is shown in Tables 1 and 2).

All data are from a single investigational device exemption center, and all surgeries were completed by 1 of 2 surgeons. Three nonrandomized Kineflex Disc training

cases were performed before beginning the randomized section of the trial; however, these data were not included in the final analysis.

There were no statistically significant differences in demographic features between the 2 groups including age, sex, body mass index, smoking status, or surgical level (Table 3). The mean preoperative VAS and ODI values were 85.4 and 63.9, respectively, for the Charité group and 83.9 and 60.1, respectively, for the Kineflex Disc group, confirming that a difference in baseline pain or dysfunction was not present and that the randomization was successful. A total of 58 of 64 patients completed at least 24 months follow-up

Table 2
Exclusion criteria

Any back or leg pain of unknown origin
 Foot drop
 Previous trauma to the study treatment level, resulting in compression or bursting
 Sufficient previous surgeries that would preclude use of an anterior approach or previous retroperitoneal surgery
 Other spinal surgery at affected level
 Previous thoracic or lumbar fusion
 Documented spondylolisthesis with >3 mm of slippage at study level
 Isthmic (spondylolytic) spondylolisthesis at study level
 Spondylitis (ie, inflammation of spine) at study level
 Documented significant spinal, foraminal, or lateral stenosis at study level
 Disc space height ≤3 mm at study level
 Documented presence of free nuclear fragment at study level
 Extensive facet arthritis or degeneration of facets at any level noted on MRI, CT, or radiography
 Scoliosis of lumbar spine with >11° of coronal deformity
 Metabolic bone disease
 Active systemic infection
 Active malignancy or history of metastatic malignancy
 Any terminal or autoimmune disease
 Any other disease, condition, or surgery that might impair healing
 Recent history of chemical or alcohol dependence
 Current or extended use of any drug known to interfere with bone or soft-tissue healing
 Known metal allergy
 Morbid obesity (body mass index >40 or >100 lb overweight)
 Transitional vertebra at level to be treated that has not clearly fused
 Pregnancy at time of enrollment (because this would contraindicate abdominal surgery)

Abbreviations: CT, computed tomography; MRI, magnetic resonance imaging.

Table 3
Study demographics

	Kineflex Disc	Charité device
No. of subjects	33	31
Male	15	13
Female	18	18
Age [mean (SD)] (y)	40.1 (10.01)	40.8 (7.68)
Body mass index [mean (SD)]	26.6 (3.39)	26.0 (3.40)
Current or former smoker	17 (51.5%)	20 (64.5%)
Surgery level L4–5	6 (18%)	4 (12.9%)
Surgery level L5-S1	27 (81.8%)	27 (87.1%)

with in the specified window (+/– two months). Three patients from the Kineflex Disc group and three from the Charité group were dropped from the study or lost before this follow-up. In the Kineflex Disc group, one patient had revision to an ALIF within 24 hours of the index operation, one patient refused follow-up, and the third had an ODI that changed from 44 to 8 and a VAS score that changed from 87 to 8 at 24 months, but these were out-of-window and thus not included. In the Charité group, one patient underwent revision to a more central location of an L5-S1 implant within 24 hours of index surgery and two patients refused follow-up. One of these was seen at 5 years' follow-up exceeding the FDA success criteria but out of window.

Primary endpoints

Several endpoints were chosen to evaluate a variety of criteria. FDA “clinical success” was assessed, defined as (1) a minimum of 25% improvement in low-back pain ODI score at 24 months compared with baseline; (2) no revisions, reoperations, or device removals; and (3) no occurrence of a major device-related adverse event.

We also compared mean improvement in VAS and ODI scores at 24 months. As cutoffs for non-inferiority, we used a VAS difference of greater than 18 points and an ODI score of greater than 10 units. No absolute consensus exists regarding what constitutes a minimal clinically important difference for the ODI and VAS scales. However, there is a good deal of literature suggesting that our chosen cutoffs are both theoretically sound and practical values that indicate that the patient has had significant relief (or worsening).^{20–24}

Non-inferiority was defined as a mean reduction in VAS score in the Kineflex Disc group at 24 months of less than 18 points below the mean reduction in the Charité group. Similarly, the Kineflex Disc device would be determined to be inferior to the Charité disc if the mean reduction in the ODI score in the Kineflex Disc group was greater than 10 units below the mean improvement in the Charité group.

Secondary endpoints

Because a major benefit of disc replacement as compared with other surgical modalities is in its long-term preservation of motion and decrease in adjacent-disc disease, we performed a mixed-effects model analysis to observe the

difference between the 2 groups over time. This model was used to understand the time course change over the follow-up period. The analysis accounts for the correlation between a series of observations taken on the same subject.

Surgical technique

The surgical technique for the Charité artificial disc implantation has been previously described in detail by Geisler et al.¹² The Kineflex artificial disc is implanted in a similar fashion, with the following differences. The Kineflex Disc comes in 3 endplate sizes. Size 1 is 27 mm × 36.5 mm, size 2 is 30 mm × 41 mm, and size 3 is 35 mm × 44 mm. The implant is inserted as a single unit with assembled disc heights of 11 mm, 12 mm, 12.75 mm, 13 mm, and 13.75 mm. In addition, there are 3 different built-in lordosis angles of 0°, 5°, and 10°. Unique to this system are 3 distractor instruments of 9 mm, 11 mm, and 13 mm in height, as well as a disc space distractor similar to the Charité distractor and a slot cutter similar to the ProDisc-L slot-cutting technique. Although the implant is 3 separate parts, it is inserted as a single unit. The core component of the implant is a constant 10.5 mm in height with a diameter of 19.25 mm and sits within the concave surfaces of the disc endplates.

Statistical analysis

Our statistical analysis was performed in accordance with the recommendations made in the CONSORT (Consolidated Standards of Reporting Trials) consensus statement for reporting of non-inferiority trials.²⁵ For our primary outcomes, a 2-sample *t* test was performed comparing the mean change in ODI and VAS scores from baseline to 2 years, as well as the percentage of each group achieving FDA-defined clinical success. Our hypothesis was that of non-inferiority; therefore the 2-sample *t* test was 1-sided. *P* values were adjusted accordingly, and the cutoff for significance (α level) was changed to .025 instead of the standard .05. Histograms were used to determine the normality of improvement scores in both the Charité and Kineflex Disc groups, for both the 2-year and 3-year analyses. Normality assumptions were adequately met for the *t* tests.

In the mixed-effects model analysis, the mean VAS pain score and mean ODI score were compared for the 2 groups over the 2-year follow-up. An interaction term between time and group was tested for each outcome to determine whether the difference between the groups changed over the follow-up period.

Results

In the Kineflex Disc group, the mean preoperative ODI and VAS scores were 60.09 and 83.89, respectively, whereas at 24 months postoperatively, the ODI and VAS scores were 22.79 and 27.09, respectively. For the Charité group, the ODI and VAS scores were 63.90 and 85.41, respectively, preoperatively but 25.50 and 30.98, respectively, at 24 months. One Charité patient underwent revision

Table 4
Mean improvement in VAS and ODI with time

	Kineflex mean improvement	Charité mean improvement	Difference (Kineflex – Charité)	95% CI	P value
VAS score					
3-y change	59.05	68.47	–9.41	–27.3 to 8.5	.17
2-y change	56.80	54.43	2.37	–12.5 to 17.3	.004
ODI					
3-y change	40.00	43.67	–3.67	–18.9 to 11.6	.2
2-y change	37.30	38.40	–1.09	–12.0 to 9.9	.054

of her L5-S1 implant to a more central position. One Kineflex patient underwent ALIF revision of her L5-S1 implant, which had subluxated anteriorly 4 mm 24 hours after the index operation. There were no other surgical or postsurgical complications, revisions, or reoperations. Every patient was discharged within 24 hours of the operation.

Table 4 shows the mean improvement for both the VAS pain and ODI scores for the Kineflex Disc and Charité arms. The mean improvement in VAS pain score (from preoperatively to 2 years postoperatively) of the Kineflex Disc group compared with the Charité group was 2.37 (95% confidence interval [CI], –12.5 to 17.3; $P = .004$). Therefore it met the assumptions of non-inferiority (difference of <18 points).

The mean improvement in ODI score from preoperatively to 24 months was –1.09 (95% CI, –12.0 to 9.9; $P = .054$). Because the CI is not constrained within the 10-point difference set as the criterion for significance, non-inferiority in terms of ODI is inconclusive.

At 24 months, 83% of patients in the Kineflex Disc group and 86% of patients in the Charité group met FDA criteria for clinical success. As described earlier, this is a composite endpoint that takes into account both relief of disability (as assessed by a greater than 25% reduction in ODI) and device safety (defined as no revisions, reoperations, removals, or major device-related adverse events). No statistically significant differences between the 2 groups existed in terms of clinical success ($P = .802$), again confirming non-inferiority.

At no point in our mixed-effects model analysis was there a significant difference between the groups in terms of mean ODI or VAS score (Figs. 1 and 2). The results from this analysis show a trend toward improved pain reduction in the Kineflex Disc group. With analysis in this fashion, the Charité arm, over time, was on average 1.02 units (95% CI, –4.10 to 6.14) higher than the Kineflex Disc group ($P = .69$). A similar trend was noted for the disability scale, where the Charité group, over time, was on average 2.43 points (95% CI, –2.33 to 7.20) higher than the Kineflex Disc group ($P = .31$).

Discussion

The purpose of a non-inferiority trial is to use a particular intervention either as a surrogate for the gold standard or in comparing two interventions where one is believed to have an important feature that is not readily amenable to a clinical trial (price, ease of use, and so on). The Kineflex Disc was designed to be introduced as a single unit, thereby aiding reliable midline posterior introduction despite a severely diseased disc space. In addition, material studies of prostheses suggest superiority of metal-on-metal devices in terms of long-term wear profiles and of reduced immunogenicity.^{26–28} However, these facts are only relevant insofar as the device can be shown to be non-inferior to an established ADR device in terms of safety and efficacy.

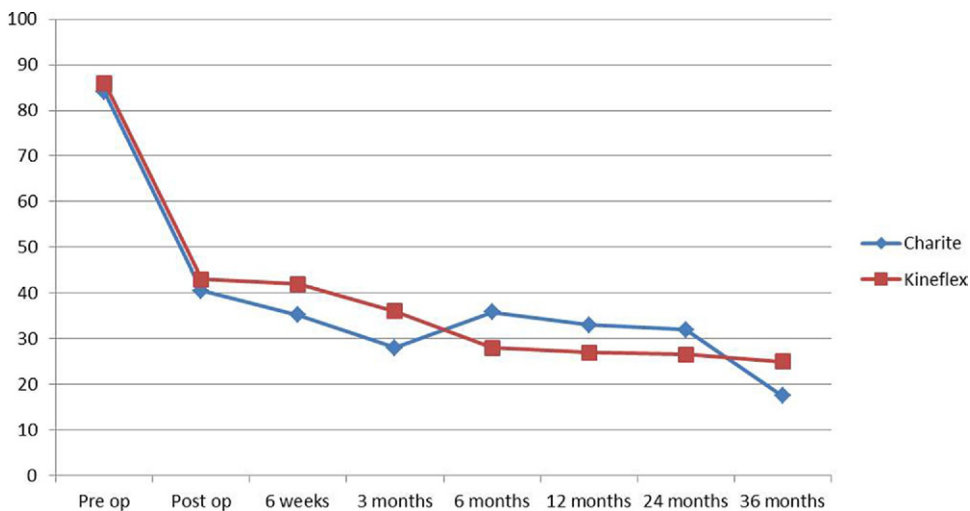


Fig. 1. Mean VAS pain scores over time (preoperatively, postoperatively, and at 6 weeks and 3, 6, 12, 24, and 36 months) for each group.

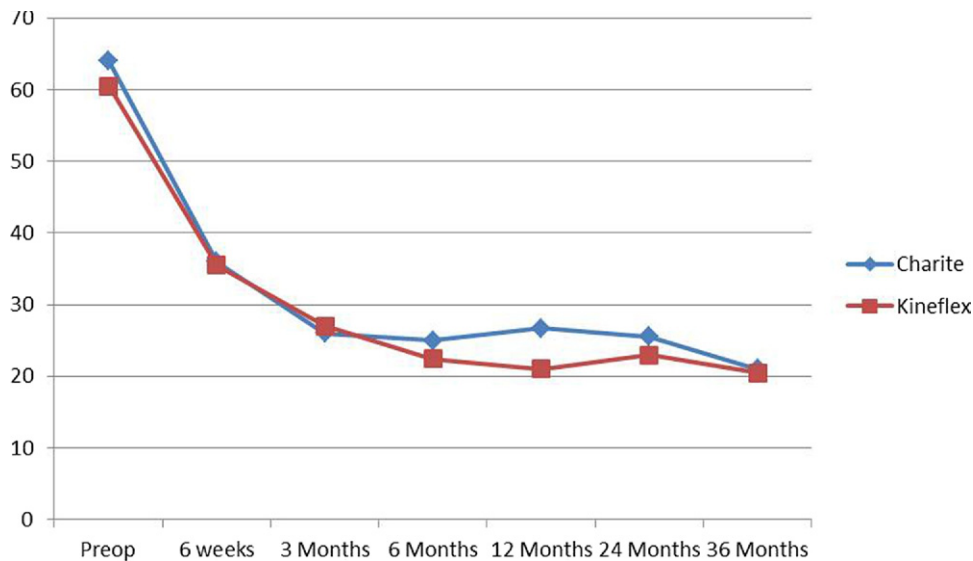


Fig. 2. Mean ODI scores over time (preoperatively and at 6 weeks and 3, 6, 12, 24, and 36 months) for each group.

In this study both groups showed statistical and clinically significant success over baseline for treatment of lower back pain caused by degenerative disc disease.

We found the Kineflex Disc to be non-inferior to the Charité device at 24 months of follow-up, defined as a mean reduction in VAS score by a difference of less than 18 points established by previous literature. The Kineflex Disc arm showed a small trend toward improvement over the Charité group (2.37 points); however, superiority cannot be claimed by this analysis.

The ODI analysis was inconclusive. This does not, however, indicate inferiority in terms of disability reduction. In fact, at 2 years, the ODI change CI includes a term nearing a minimum clinically significant improvement.

In addition, the mixed-effects model analysis confirms Kineflex Disc non-inferiority over time and shows a small trend toward superiority in pain and disability scores, although these conclusions cannot be definitively made.

Limitations

Our study is potentially limited by the problem of “bio-creep” inherent in any multigenerational non-inferiority trial.²⁹ With this problem, each subsequent iteration of non-inferiority trial increases the size of error bars around the true results, and it becomes increasingly difficult to state with certainty that the newest intervention is, in fact, non-inferior to the original (in this case ALIF with BAK cages). That said, because this procedure is no longer the standard of care, such a trial is not possible despite it being preferable in terms of statistical modeling.

In addition, a conclusion that one might be tempted to draw from this study (ie, non-inferiority of the Kineflex Disc device to spinal fusion) is limited by the fact that original trials comparing the Charité device with fusion used ALIF with BAK cages, a procedure no longer considered to be the gold standard for fusions. However,

because this was the standard of care at the time of the original Charité trials, this control arm was considered reasonable.

Conclusions

The Kineflex Disc ADR device is non-inferior to the Charité ADR device at 24 months’ follow-up in terms of pain reduction and FDA-defined clinical success—a composite term that includes safety measures. The ease and reliability of implantation of the Kineflex Disc ADR device and improved wear profiles make it a viable alternative technology. Additional studies are needed to further characterize the long-term risks and benefits of these devices beyond 2-year follow-up, and a well-designed head-to-head trial is needed to establish whether any particular ADR is superior with respect to hard clinical outcomes. Every patient (except the 2 undergoing reoperations) was discharged from the hospital within 24 hours of surgery. These cases are scheduled as outpatient procedures.

References

1. Yajun W, Yue Z, Xiuxin H, Cui C. A meta-analysis of artificial total disc replacement versus fusion for lumbar degenerative disc disease. *Eur Spine J* 2010;19:1250–61.
2. Gamradt SC, Wang JC. Lumbar disc arthroplasty. *Spine J* 2005;5:95–103.
3. McAfee PC, Yu W, Wang JC. Controversies in surgical treatment for LBP: fusion or disc replacement. *SpineLine* 2004:14–20.
4. Berg S, Tullberg T, Branth B, Olerud C, Tropp H. Total disc replacement compared to lumbar fusion: a randomised controlled trial with 2-year follow-up. *Eur Spine J* 2009;18:1512–9.
5. Blumenthal S, McAfee PC, Guyer RD, et al. A prospective, randomized, multicenter Food and Drug Administration investigational device exemptions study of lumbar total disc replacement with the CHARITE artificial disc versus lumbar fusion: part I: evaluation of clinical outcomes. *Spine* 2005;30:1565–75.

6. Delamarter RB, Bae HW, Pradhan BB. Clinical results of ProDisc-II lumbar total disc replacement: report from the United States clinical trial. *Orthop Clin North Am* 2005;36:301–13.
7. Sasso RC, Foulk DM, Hahn M. Prospective, randomized trial of metal-on-metal artificial lumbar disc replacement: initial results for treatment of discogenic pain. *Spine* 2008;33:123–31.
8. Harrop JS, Youssef JA, Maltenfort M, et al. Lumbar adjacent segment degeneration and disease after arthrodesis and total disc arthroplasty. *Spine* 2008;33:1701–7.
9. Gillet P. The fate of the adjacent motion segments after lumbar fusion. *J Spinal Disord Tech* 2003;16:338–45.
10. Kumar MN, Jacquot F, Hall H. Long-term follow-up of functional outcomes and radiographic changes at adjacent levels following lumbar spine fusion for degenerative disc disease. *Eur Spine J* 2001;10:309–13.
11. Park P, Garton HJ, Gala VC, Hoff JT, McGillicuddy JE. Adjacent segment disease after lumbar or lumbosacral fusion: review of the literature. *Spine* 2004;29:1938–44.
12. Geisler FH, Blumenthal SL, Guyer RD, et al. Neurological complications of lumbar artificial disc replacement and comparison of clinical results with those related to lumbar arthrodesis in the literature: results of a multicenter, prospective, randomized investigational device exemption study of Charité intervertebral disc. Invited submission from the Joint Section Meeting on Disorders of the Spine and Peripheral Nerves, March 2004. *J Neurosurg Spine* 2004;1:143–54.
13. McAfee PC, Cunningham B, Holsapple G, et al. A prospective, randomized, multicenter Food and Drug Administration investigational device exemption study of lumbar total disc replacement with the CHARITE artificial disc versus lumbar fusion: part II: evaluation of radiographic outcomes and correlation of surgical technique accuracy with clinical outcomes. *Spine* 2005;30:1576–83.
14. Griffith SL, Shelokov AP, Büttner-Janzen K, LeMaire JP, Zeegers WS. A multicenter retrospective study of the clinical results of the LINK SB Charité intervertebral prosthesis. The initial European experience. *Spine* 1994;19:1842–9.
15. Lemaire JP, Carrier H, Soriali el-H, Skalli W, Lavaste F. Clinical and radiological outcomes with the Charité artificial disc: a 10-year minimum follow-up. *J Spinal Disord Tech* 2005;16:353–9.
16. Guyer RD, McAfee PC, Banco RJ, et al. Prospective, randomized, multicenter Food and Drug Administration investigational device exemption study of lumbar total disc replacement with the CHARITE artificial disc versus lumbar fusion: five-year follow-up. *Spine J* 2009;9:374–86.
17. Khatri M, Norris H, Ross ERS. Results of 226 Charite III arthroplasties done in 160 patients with an average follow-up of 79 month. Presented at the Annual Meeting of the International Society for the Study of the Lumbar Spine, New York, NY, May 2005.
18. Geisler FH. Surgical treatment of discogenic low-back pain: lumbar arthroplasty results in superior pain reduction and disability level improvement compared with lumbar fusion. *SAS J* 2007;1:12–9.
19. Hanle UR, Weinberg IR, Sliwa K, Sweet BM, de Villiers M. Kineflex (Centurion) lumbar disc prosthesis: insertion technique and 2-year clinical results in 100 patients. *SAS J* 2007;1:28–35.
20. Hägg O, Fritzell P, Nordwall A, Swedish Lumbar Spine Study Group. The clinical importance of changes in outcome scores after treatment for chronic low back pain. *Eur Spine J* 2003;12:12–21.
21. Maughan EF, Lewis JS. Outcome measures in chronic low back pain. *Eur Spine J* 2010;19:1484–94.
22. Copay AG, Glassman SD, Subach BR, et al. Minimum clinically important difference in lumbar spine surgery patients: a choice of methods using the Oswestry Disability Index, Medical Outcomes questionnaire Short Form 36, and Pain Scales. *Spine J* 2009;8:968–74.
23. Copay AG, Martin MM, Subach BR, et al. Assessment of spine surgery outcomes: inconsistency of change amongst outcome measurements. *Spine J* 2010;10:291–6.
24. Jensen MP, Chen C, Brugger AM. Interpretation of visual analog scale ratings and change scores: a reanalysis of two clinical trials of post-operative pain. *J Pain* 2003;4:407–14.
25. Piaggio G, Elbourne DR, Altman DG, Pocock SJ, Evans SJ, CONSORT Group. Reporting of noninferiority and equivalence randomized trials: an extension of the CONSORT statement. *JAMA* 2006;295:1152–60.
26. van Ooij A, Kurtz SM, Stessels F, Noten H, van Rhijn L. Polyethylene wear debris and long-term clinical failure of the Charité disc prosthesis: a study of 4 patients. *Spine* 2007;32:223–9.
27. Baumann B, Rader CP, Seufert J, et al. Effects of polyethylene and TiAlV wear particles on expression of RANK, RANKL and OPG mRNA. *Acta Orthop Scand* 2004;75:295–302.
28. Kurtz SM, van Ooij A, Ross R, et al. Polyethylene wear and rim fracture in total disc arthroplasty. *Spine J* 2007;7:12–21.
29. Mirza S. The non-inferiority principle: design considerations relevant to noninferiority trials of artificial disc replacements. *SpineLine* 2006:17–20.